


THE INNOVATIVE DELTA PURPLE PICC RANGE



The catheter is a next-generation Polyurethane purple color for an easy recognition as High Pressure Resistant 300 psi.

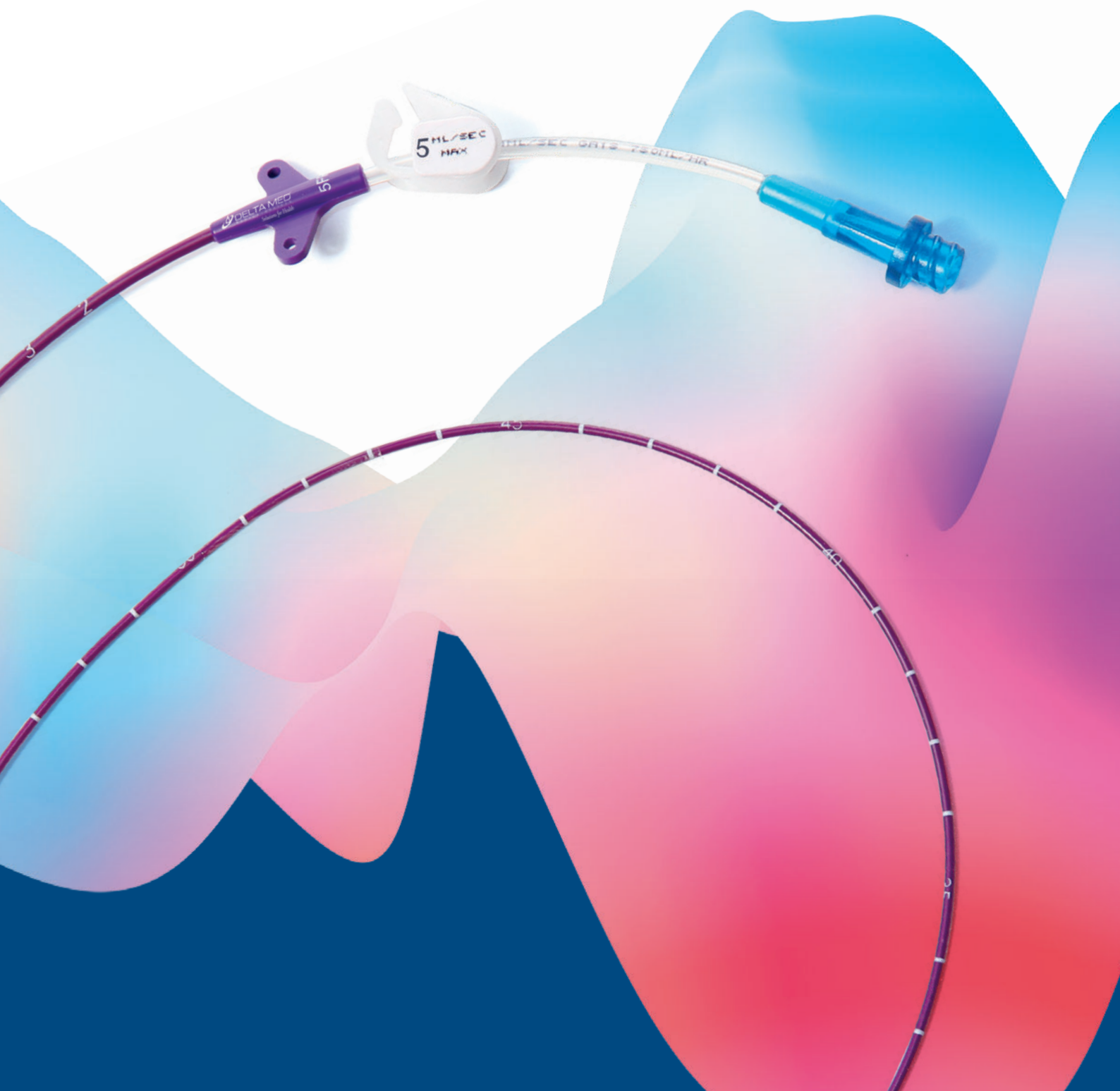

OUR SOLUTION


WEALTH OF ACCESSORIES

- Wide range: Single lumen 3/4/5 Fr; Double Lumen 4/5/6 Fr; Triple Lumen 6 Fr
- Reverse Taper
- Resistant to high pressure contrast media 300 psi
- Depth marking that can be cut distally
- Open tip device
- Equipped with T-connector support with side valve to reduce any malfunctions compared to traditional clamps during pre-washing phase.
- Presence of a rigid centimeter stylet with soft tip to facilitate the movement of the catheter during the positioning phases.
- Presence of a micro-introducer kit with echogenic safety needle 21x7cm and peel-away 7cm introducer for an easier placement without vessel trauma.
- Presence of a NEUTRON[®] quality needle less neutral connector for the washing phases after placement.
- Kit equipped with the BEDAL 2, suture less universal PICC fixing system and 10x12cm transparent polyurethane dressing for the protection of the insertion site and the fixing system after the placement phases

The catheter is indicated for the infusion of drugs for chemotherapy, antibiotic therapy, perfusions, parenteral nutrition, infusion of blood products and for intravenous blood sampling.

THE DELTA PURPLE PICC RANGE IS AN EXCELLENT SOLUTION TO INFUSE CONTRAST MEDIA WITH A MAXIMUM PRESSURE OF 300 PSI AND A FLOW RATE EQUAL TO 5 ML / SEC.

OUR BASIC KIT

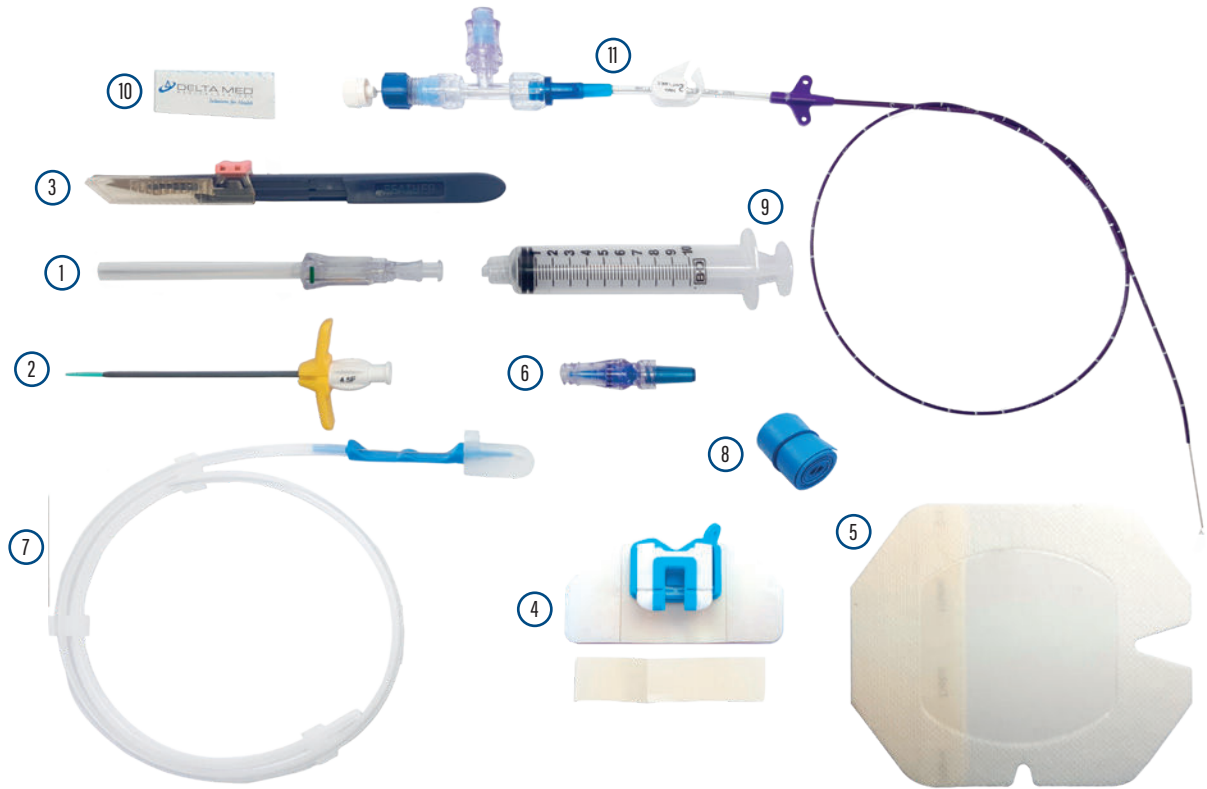

- | | | | |
|---|------------------------------------|---|-----------------------|
| ① ECHOGENIC SAFETY NEEDLE 21GX7CM | ④ FIXATION DRESSING SYSTEM | ⑦ STAINLESS STEEL GUIDEWIRE TIP 0,018"X50CM | ⑩ CENTIMETER METER |
| ② TEAR AWAY INTRODUCER 7 CM WITH INTRODUCER | ⑤ MEDICATION FOR THE STABILIZATION | ⑧ TOURNIQUET | ⑪ DELTA PICC CATHETER |
| ③ SAFETY SCALPEL | ⑥ NEUTRON NEEDLELESS CONNECTOR | ⑨ 10ML LL SYRINGE | |

OUR PLUS

MAIN FEATURES

- Extremely smooth surface = lower adhesion to the vessel wall
- Same softness as silicone elastomers
- Higher mechanic resistance than silicone
- High tensile and chemical resistance
- Greater internal diameter with the same external diameter
- Excellent biocompatibility
- Excellent thermo-sensitivity
- Demp free

LATEX FREE

ULTRASOUND SYSTEM AND TIP LOCATION

System composed of medical Tablet, Ultrasound and ECG including a wide range of Linear, Convex and Endocavitary probes from 2 to 15 MHz. Available upon request and customizable.



PICC KIT

1 | KIT CENTRAL VENOUS CATHETER IN POWER INJECTABLE POLYURETHANE WITH "PICC" PERIPHERAL INSERTION

CODE	DESCRIPTION	Q.TY	KIT COMPOSITION	
MD 01 02288/AMC	PRSLPICC-0355-T 1L 3 Fr. x 55 cm. + int. 3,5 Fr. x 7 cm	1	1 ECHOGENIC SAFETY NEEDLE 21GX7CM	9 10ML LL SYRINGE
MD 01 02297/AMC	PRSLPICC-0355-T 1L 4 Fr. x 55 cm. + int. 4,5 Fr. x 7 cm	1	2 TEAR AWAY INTRODUCER 7 CM WITH INTRODUCER	10 CENTIMETER METER
MD 01 02298/AMC	PRSLPICC-0355-T 1L 5 Fr. x 55 cm. + int. 5,5 Fr. x 7 cm	1	3 SAFETY SCALPEL	11 DELTA PICC CATHETER
MD 01 02289/AMC	PRSLPICC-0355-T 2L 4 Fr. x 55 cm. + int. 4,5 Fr. x 7 cm	2	4 FIXATION DRESSING SYSTEM	
MD 01 02295/AMC	PRSLPICC-0355-T 2L 5 Fr. x 55 cm. + int. 5,5 Fr. x 7 cm	2	5 MEDICATION FOR THE STABILIZATION	
MD 01 02296/AMC	PRSLPICC-0355-T 2L 6 Fr. x 55 cm. + int. 6,5 Fr. x 7 cm	2	6 NEUTRON NEEDLELESS CONNECTOR	
MD 01 02299/AMC	PRSLPICC-0355-T 3L 6 Fr. x 55 cm. + int. 6,5 Fr. x 7 cm	3	7 STAINLESS STEEL GUIDEWIRE TIP 0,018"X50CM	
			8 TOURNIQUET	

PICC PLUS KIT FOR ECG TECHNIQUE

1 | KIT CENTRAL VENOUS CATHETER IN POWER INJECTABLE POLYURETHANE WITH "PICC" PERIPHERAL INSERTION

CODE	DESCRIPTION	Q.TY	KIT COMPOSITION	
MD 01 02288/ASI	PRSLPICC-0355-T 1L 3 Fr. x 55 cm. + int. 3,5 Fr. x 7 cm	1	1 ECHOGENIC SAFETY NEEDLE 21GX7CM	9 10ML LL SYRINGE
MD 01 02297/ASI	PRSLPICC-0355-T 1L 4 Fr. x 55 cm. + int. 4,5 Fr. x 7 cm	1	2 TEAR AWAY INTRODUCER 7 CM WITH INTRODUCER	10 CENTIMETER METER
MD 01 02298/ASI	PRSLPICC-0355-T 1L 5 Fr. x 55 cm. + int. 5,5 Fr. x 7 cm	1	3 SAFETY SCALPEL	11 DELTA PICC CATHETER
MD 01 02289/ASI	PRSLPICC-0355-T 2L 4 Fr. x 55 cm. + int. 4,5 Fr. x 7 cm	2	4 FIXATION DRESSING SYSTEM	
MD 01 02295/ASI	PRSLPICC-0355-T 2L 5 Fr. x 55 cm. + int. 5,5 Fr. x 7 cm	2	5 MEDICATION FOR THE STABILIZATION	
MD 01 02296/ASI	PRSLPICC-0355-T 2L 6 Fr. x 55 cm. + int. 6,5 Fr. x 7 cm	2	6 NEUTRON NEEDLELESS CONNECTOR	
MD 01 02299/ASI	PRSLPICC-0355-T 3L 6 Fr. x 55 cm. + int. 6,5 Fr. x 7 cm	3	7 STAINLESS STEEL GUIDEWIRE TIP 0,018"X50CM	
			8 TOURNIQUET	

PICC PROCEDURAL KIT FOR ECG TECHNIQUE

1 | KIT CENTRAL VENOUS CATHETER IN POWER INJECTABLE POLYURETHANE WITH "PICC" PERIPHERAL INSERTION

CODE	DESCRIPTION	Q.TY	KIT COMPOSITION	
MD 01 02266/ASA	PRSLPICC-0355-T 1L 3 Fr. x 55 cm. + int. 3,5 Fr. x 7 cm	1	1 ECHOGENIC SAFETY NEEDLE 21GX7CM	9 10ML LL SYRINGE
MD 01 02260/ASA	PRSLPICC-0355-T 1L 4 Fr. x 55 cm. + int. 4,5 Fr. x 7 cm	1	2 TEAR AWAY INTRODUCER 7 CM WITH INTRODUCER	10 CENTIMETER METER
MD 01 02261/ASA	PRSLPICC-0355-T 1L 5 Fr. x 55 cm. + int. 5,5 Fr. x 7 cm	1	3 SAFETY SCALPEL	11 DELTA PICC CATHETER
MD 01 02265/ASA	PRSLPICC-0355-T 2L 4 Fr. x 55 cm. + int. 4,5 Fr. x 7 cm	2	4 FIXATION DRESSING SYSTEM	
MD 01 02262/ASA	PRSLPICC-0355-T 2L 5 Fr. x 55 cm. + int. 5,5 Fr. x 7 cm	2	5 MEDICATION FOR THE STABILIZATION	
MD 01 02263/ASA	PRSLPICC-0355-T 2L 6 Fr. x 55 cm. + int. 6,5 Fr. x 7 cm	2	6 NEUTRON NEEDLELESS CONNECTOR	
MD 01 02264/ASA	PRSLPICC-0355-T 3L 6 Fr. x 55 cm. + int. 6,5 Fr. x 7 cm	3	7 STAINLESS STEEL GUIDEWIRE TIP 0,018"X50CM	
			8 TOURNIQUET	

1 pc. 3 ml LI Syringe - 1 pc. 5 ml LL Syringe
 1 pc. 10 ml LL Syringe - 1 pc. Safety hypodermic needle 25x25mm - 1 pc. Probe cover sheath 15x122cm and single-dose gel - 1 pc. IV dressing 10x12 cm with reinforced edge - 1 pc. Guillotine for catheter cutting - 10 pcs. Cotton gauzes 10x10 cm 8 layers - 1 pc. >Blue rectangula tray 40x26x7cm - 1 pc. Drape in SMS 180x330cm with right anmd left arm holes - 1 pc. Double-coupled drape 75x90xm - 1 pc. Wrap-around surgical gown in SMS size L with 2 wipes - 1 pc. Cap with elastic for hair - 1 pc. Surgical mask with elastics - 1 Pair of powder-free synthetic surgical gloves size 7.5 ULTRAPRENE SYNTHETIC - 1 pc. Centimeter Meter "tape of paper" - 1 pc. Tourniquet - enveloping hydro alcohol repellent SMS towel.