

Distributed by

# Jackson-Pratt®

# The name you trust for wound drainage products



Wound Drainage Products

# Jackson-Pratt® Wound Drainage Products

For over 40 years, Jackson-Pratt® has been a leading name in the wound drainage industry. Delivering high-quality, reliable products has allowed the Cardinal Health Jackson-Pratt® brand to become a trusted solution worldwide.

When you feel confident about the products you use, you can focus on what matters most — providing the best care possible to your patients.

## **Jackson-Pratt®**

The name you trust for wound drainage products.

3-4 **Reservoirs and Accessories**

5-6 **Silicone and PVC Wound Drains**

7 **Hemaduct® Wound Drains**






100cc and 400cc bulb reservoirs

## Jackson-Pratt® 100cc and 400cc Silicone Bulb Reservoirs and Reservoir Kits

For procedure needs, Jackson-Pratt® 100cc and 400cc Reservoirs are available separately for use with all Jackson-Pratt® Wound Drains or in convenient kits packaged with Jackson-Pratt® Flat or Round Wound Drains.

### Jackson-Pratt® Reservoirs feature:

- Silicone walls to provide easy activation of suction
- Internal anti-reflux valve to help prevent backflow of fluid to the patient
- Clear reservoir and graduations that make it easy to identify and measure exudate
- Not made with natural rubber latex
- 100cc bulb that features a Luer lock emptying port which connects directly to either a Luer lock syringe for collection of a culture sample or to an exudate disposal bag to provide a closed system for safe disposal of exudate
- 400cc bulb that features dual-inlet ports to allow connection of two drains without a y-connector
- 400cc bulb with a threaded outlet port that connects directly to an exudate disposal bag for safe disposal of exudate

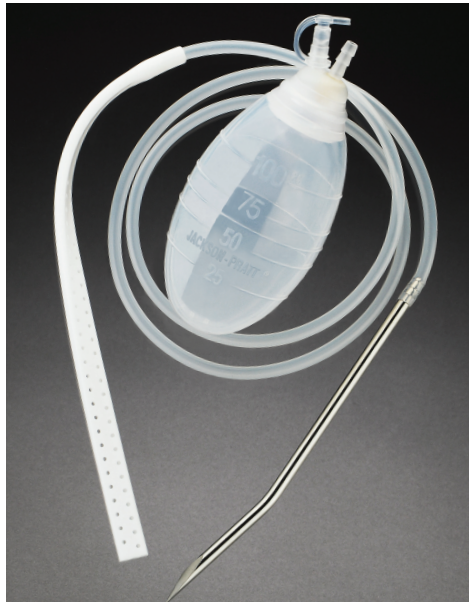

100cc reservoir kit with flat silicone drain with trocar



Bifurcated drain extension, drain to non-Jackson-Pratt® Reservoir Connector and wall suction connector with tubing to connect to 100cc reservoir


## Jackson-Pratt® 100cc and 400cc Silicone Bulb Reservoirs and Reservoir Kits (continued)

Cat. No.	Description	Perforation
<b>100cc silicone bulb reservoir (30 ea/cs)</b>		
SU130-1305D	100cc reservoir only	n/a
<b>100cc reservoir kits with flat silicone drains without trocar (10/bx)</b>		
SU130-1348	7mm x 20cm	¾
SU130-1349	10mm x 20cm	¾
SU130-1360	7mm x 20cm	full
SU130-1361	10mm x 20cm	full
<b>100cc reservoir kits with flat silicone drains with trocars (10/bx)</b>		
SU130-1358	7mm x 20cm	¾
SU130-1359	10mm x 20cm	¾
SU130-1370	7mm x 20cm	full
SU130-1371	10mm x 20cm	full

Cat. No.	Description	Perforation
<b>400cc silicone bulb reservoir (10/bx)</b>		
SU130-1000	400cc reservoir only	n/a
<b>400cc silicone bulb reservoir (10EA/CS)</b>		
SU1301000D	400cc reservoir only	n/a
<b>400cc reservoir kits with round silicone drains with trocar (10/bx)</b>		
SU130-1080	10 Fr, 1/8"	center
SU130-1060	15 Fr, 3/16"	end
SU130-1010	19 Fr, ¼"	end

## Jackson-Pratt® Accessories

Cat. No.	Description	Size	Qty.
SU130-1315	Bifurcated drain extension, accommodates dual drainage, as required	5mm O.D.	10/bx
SU130-1340	Drain to non-Jackson-Pratt® Reservoir Connector	3mm to 6mm	10/bx
SU130-1342	Wall suction connector, with 2" length, 6mm (¼") tubing to connect to 100cc reservoir	3mm to 10mm	10/bx

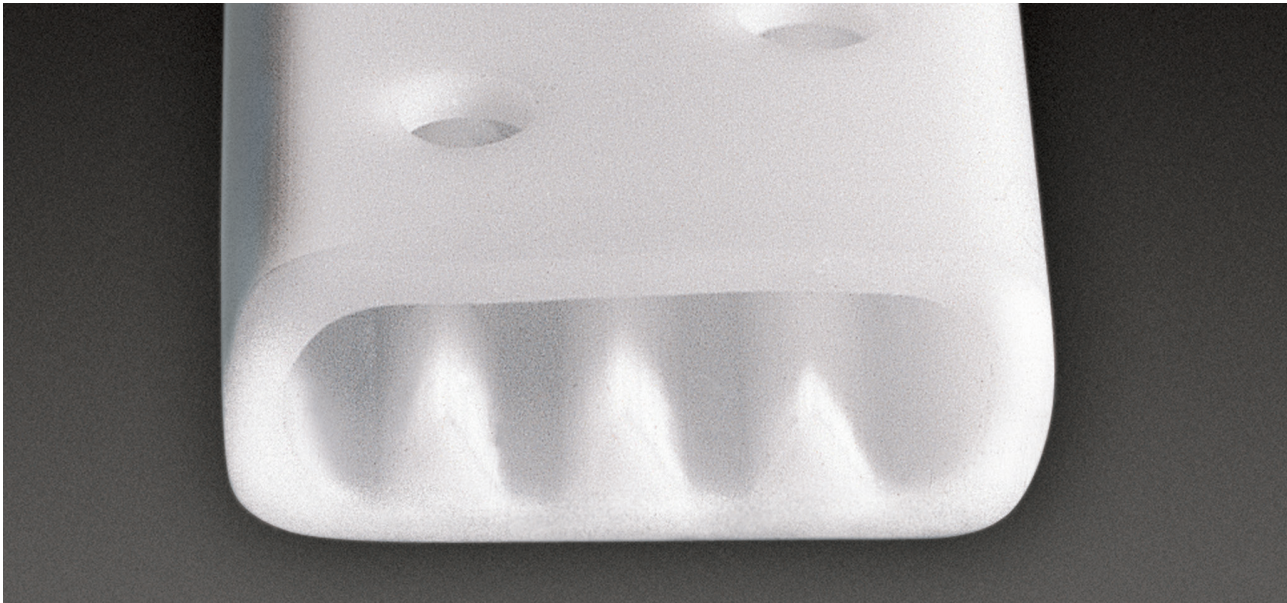

Jackson-Pratt® Silicone Flat Drain

## Jackson-Pratt® Silicone Flat Drains

**Jackson-Pratt® Flat Drains are made of a silicone elastomer for softness and flexibility, featuring:**

- Injection molded for strength and drain integrity
- Low-profile, hubless design to help promote tissue plane approximation
- Inner lumen ribbing to help prevent drain collapse and clogging and help preserve drain patency
- Silicone that is barium impregnated for x-ray detection of drain

Cat. No.	Drain size	Perforation
<b>Flat silicone drains without trocar (10/bx)</b>		
SU130-1308	7mm x 20cm	¾
SU130-1309	10mm x 20cm	¾
SU130-1310	7mm x 20cm	full
SU130-1311	10mm x 20cm	full
<b>Flat silicone drains with 15 Fr trocar (10/bx)</b>		
SU130-1410	7mm x 20cm	full
SU130-1411	10mm x 20cm	full

# Silicone and PVC Wound Drains



Silicone round drain with trocar



Silicone t-tube drain

## Jackson-Pratt® Silicone and PVC Round Drains

Jackson-Pratt® Round Drains are available in both silicone and PVC per physician preference and feature a radiopaque stripe running the entire length of the drain tubing for x-ray detection.

Cat. No.	Drain size	Perforation
<b>Silicone round drains without trocar (10/bx)</b>		
SU130-1320	7 Fr (3/32")	end
SU130-1321	10 Fr (1/8")	end
SU130-1323	15 Fr (3/16")	end
SU130-1325	19 Fr (1/4")	end
SU130-1327	19 Fr (1/4")	center
SU130-1328	10 Fr (1/8")	center
SU130-1329	15 Fr (3/16")	center

<b>Silicone round drains with trocar (10/bx)</b>		
SU130-0320	7 Fr (3/32")	end
SU130-0321	10 Fr (1/8")	end
SU130-0323	15 Fr (3/16")	end
SU130-0325	19 Fr (1/4")	end

<b>PVC round drains with trocar (10/bx)</b>		
SU130-0521	10 Fr (1/8")	center
SU130-0523	15 Fr (3/16")	center
SU130-0525	19 Fr (1/4")	center


## Jackson-Pratt® Silicone T-tube Drain

The Jackson-Pratt® Silicone T-tube Drain is indicated for hysterectomy and cholecystectomy and is soft and pliable.

Cat. No.	Description	Drain Size	Qty.
SU130-1334	Silicone t-tube drain	19 Fr (1/4")	10/bx



Hemaduct® Flat Drain

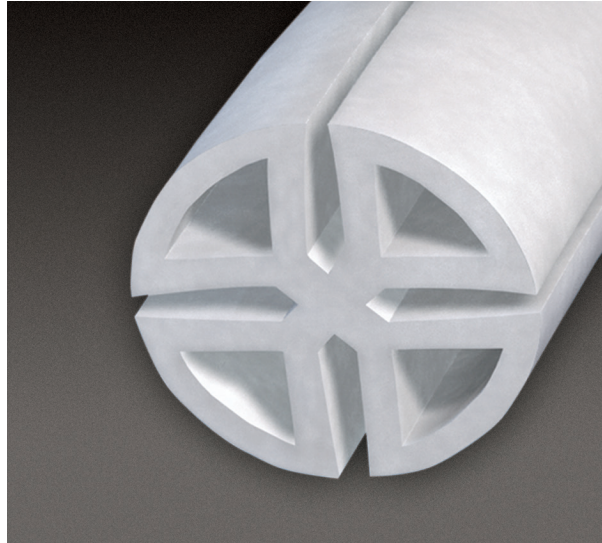

Hemaduct® Round Drain

The Jackson-Pratt® brand Hemaduct® Wound Drain is uniquely designed with a system of multiple ducts (channels) and lumens that are interconnected through a series of internal portals.

**This advanced design offers the following advantages:**

- Provides suction and fluid flow over the entire drain
- Provides alternate pathways around obstructions
- Hemaduct® Drains are made of soft silicone and are available in round and flat styles, with or without a trocar

Cat. No.	Drain Size*	Ducts	Trocar Size
----------	-------------	-------	-------------

**Hemaduct® flat drains without trocar (10/bx)**

JP-HUF070	7mm x 20cm	Full (20cm)	
JP-HUF100	10mm x 20cm	Full (20cm)	

**Hemaduct® flat drains with trocar (10/bx)**

JP-HUF071	7mm x 20cm	Full (20cm)	15 Fr x 15cm
JP-HUF104	10mm x 20cm	¾ (15cm)	15 Fr x 15cm
JP-HUF101	10mm x 20cm	Full (20cm)	15 Fr x 15cm

**Hemaduct® ¾ ducted round drains without trocar (10/bx)**

JP-HUR100	10 Fr (1/8") x 40cm	30cm	
JP-HUR150	15 Fr (3/16") x 40cm	30cm	
JP-HUR190	19 Fr (¼") x 40cm	30cm	

**Hemaduct® ¾ ducted round drains with trocar (10/bx)**

JP-HUR101	10 Fr (1/8") x 40cm	30cm	10 Fr x 15cm
JP-HUR151	15 Fr (3/16") x 40cm	30cm	15 Fr x 15cm
JP-HUR195	19 Fr (¼") x 40cm	30cm	15 Fr x 15cm

\*Length from hub to distal end

Cat. No.	Drain Size*	Ducts	Trocar Size
----------	-------------	-------	-------------

**Hemaduct® fully ducted round drains without trocar (10/bx)**

JP-HUR860	10 Fr (1/8") x 30cm	30cm	
JP-HUR880	15 Fr (3/16") x 30cm	30cm	
JP-HUR900	19 Fr (¼") x 30cm	30cm	

**Hemaduct® fully ducted round drains with trocar (10/bx)**

JP-HUR870	10 Fr (1/8") x 30cm	30cm	10 Fr x 15cm
JP-HUR890	15 Fr (3/16") x 30cm	30cm	15 Fr x 15cm
JP-HUR910	19 Fr (¼") x 30cm	30cm	15 Fr x 15cm

# Hemaduct® Wound Drains

This literature piece is for distribution through ISKUS HEALTH in the United Kingdom and Ireland only.

© 2015 Cardinal Health. All Rights Reserved. CARDINAL HEALTH, the Cardinal Health LOGO, ESSENTIAL TO CARE, HEMADUCT and JACKSON-PRATT are trademarks or registered trademarks of Cardinal Health. All other marks are the property of their respective owners. Lit. No. 2PERI15-370597 (08/2015) PRC (08/2015)



Distributed by



[www.iskushealth.com](http://www.iskushealth.com)

#### Customer Services and Distribution Contact Details

##### Ireland

Iskus Health Ltd  
21 Fonthill Business Park  
Fonthill Road  
Dublin 22

**Tel:** +353 1 404 8383

**Fax:** +353 1 404 8381

##### From Northern Ireland

**Tel:** 0845 850 8600

**Fax:** 0845 850 8601

**Email:** [info@iskushealth.com](mailto:info@iskushealth.com)

##### United Kingdom

Iskus Health UK Ltd.  
B1-B2 Bond Close  
Kingsland Business Park  
Basingstoke  
RG24 8PZ

**Tel:** +44 (0) 1256365458

**Fax:** +44 (0) 1256365459

**Email:** [uksales@iskushealth.com](mailto:uksales@iskushealth.com)